



4 Main Features for New OF4 Series

Forced_Convection_Oven



4 Main Features for New OF4 Series

01

5-Inch TFT touch screen controller 🔍



OF3 Series
(VFD Display Controller)



NEW
OF4 Series

- ✓ Excellent visibility and intuitive control with a 5-inch touch screen.
- ✓ Three product lines classified by different controller types(OF4-P/V/S)
- ✓ Improved functionality such as real-time graphs and ramping rate setting

4 Main Features for New OF4 Series

02

3-Step Fan Speed Adjustment




Opti-flow™ function


Even if Forced Convection Oven, it allows prevent floating dust particles in the chamber like a Natural convection oven through **flexibly adjustable fan speed** according to the specimen type and test conditions.

Drying of samples that can cause dust

Opti-flow™ level 1

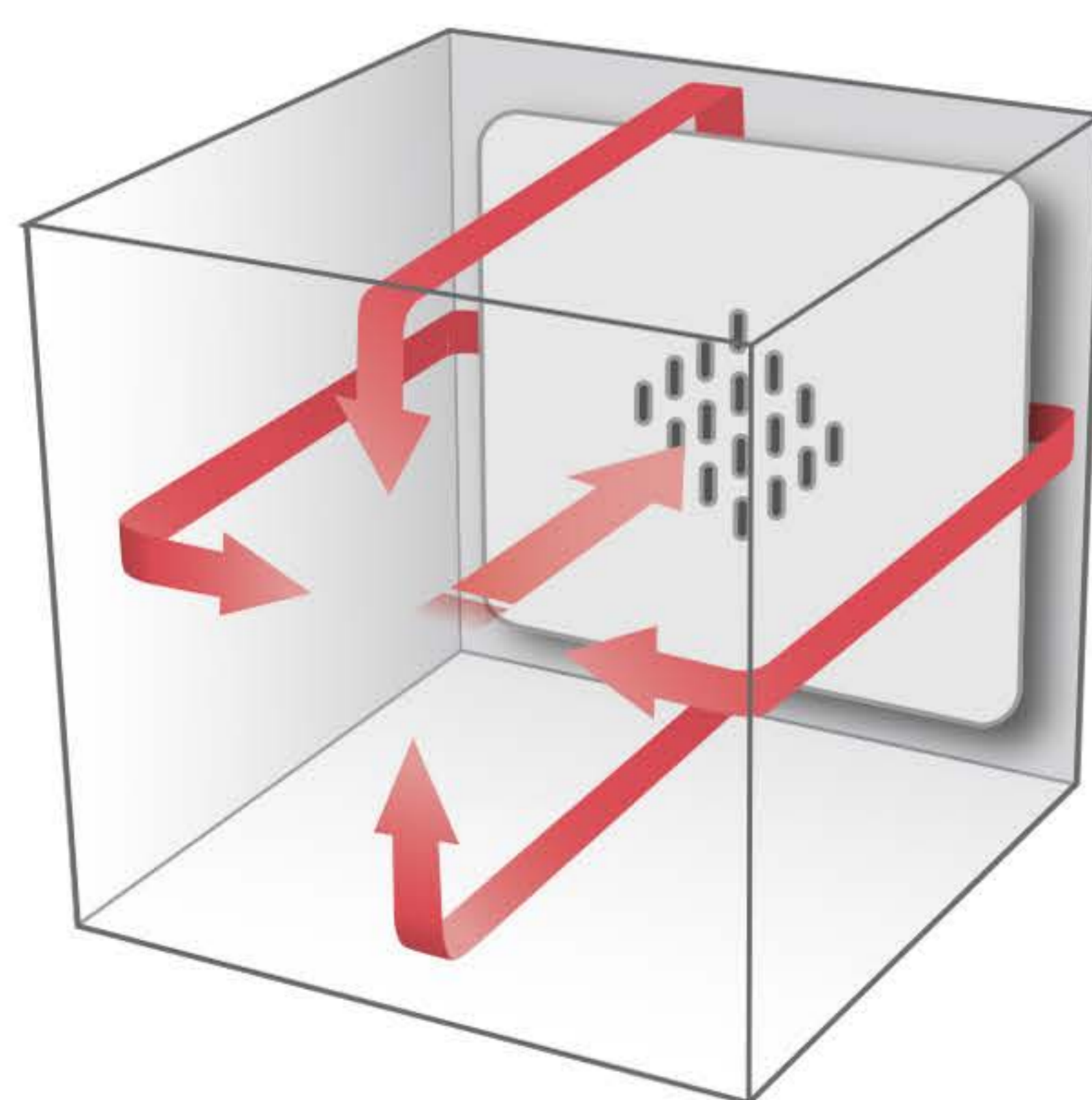


Opti-flow™ level 3

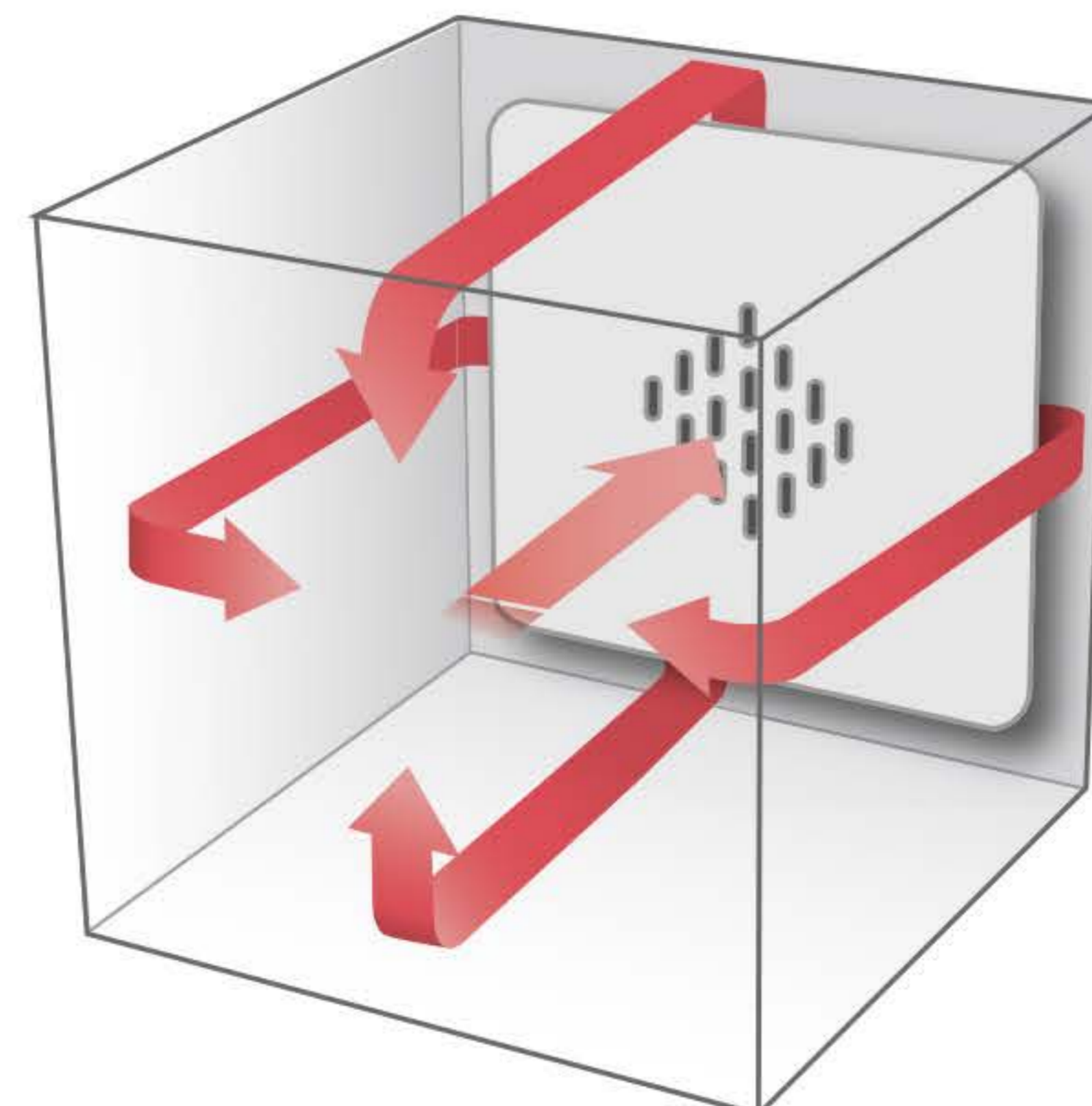

Drying at the fastest rate if dust isn't generated

Opti-flow™ level 5

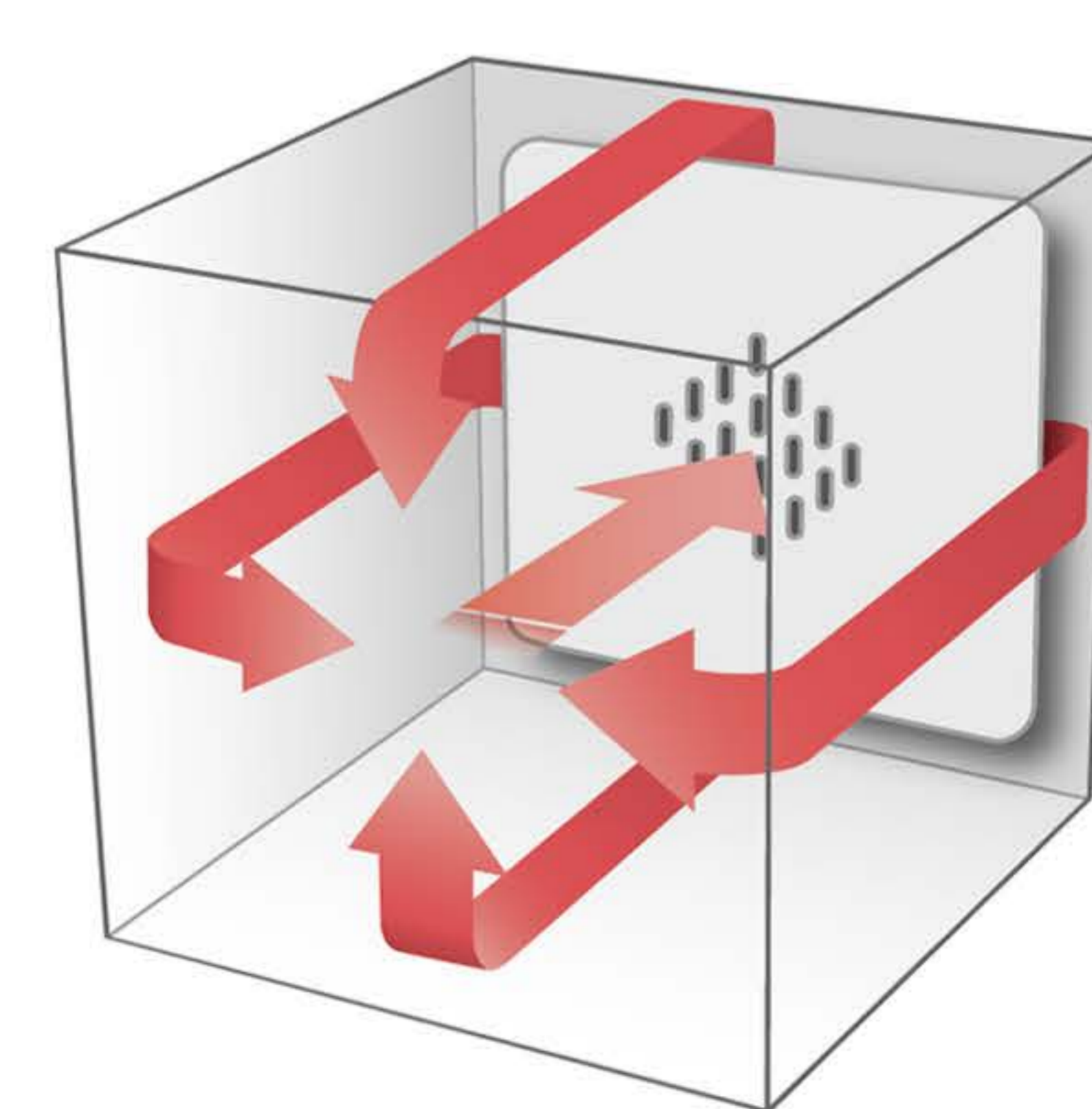


4 Main Features for New OF4 Series

03

Convenient Data Management



OF-4 V/P Model

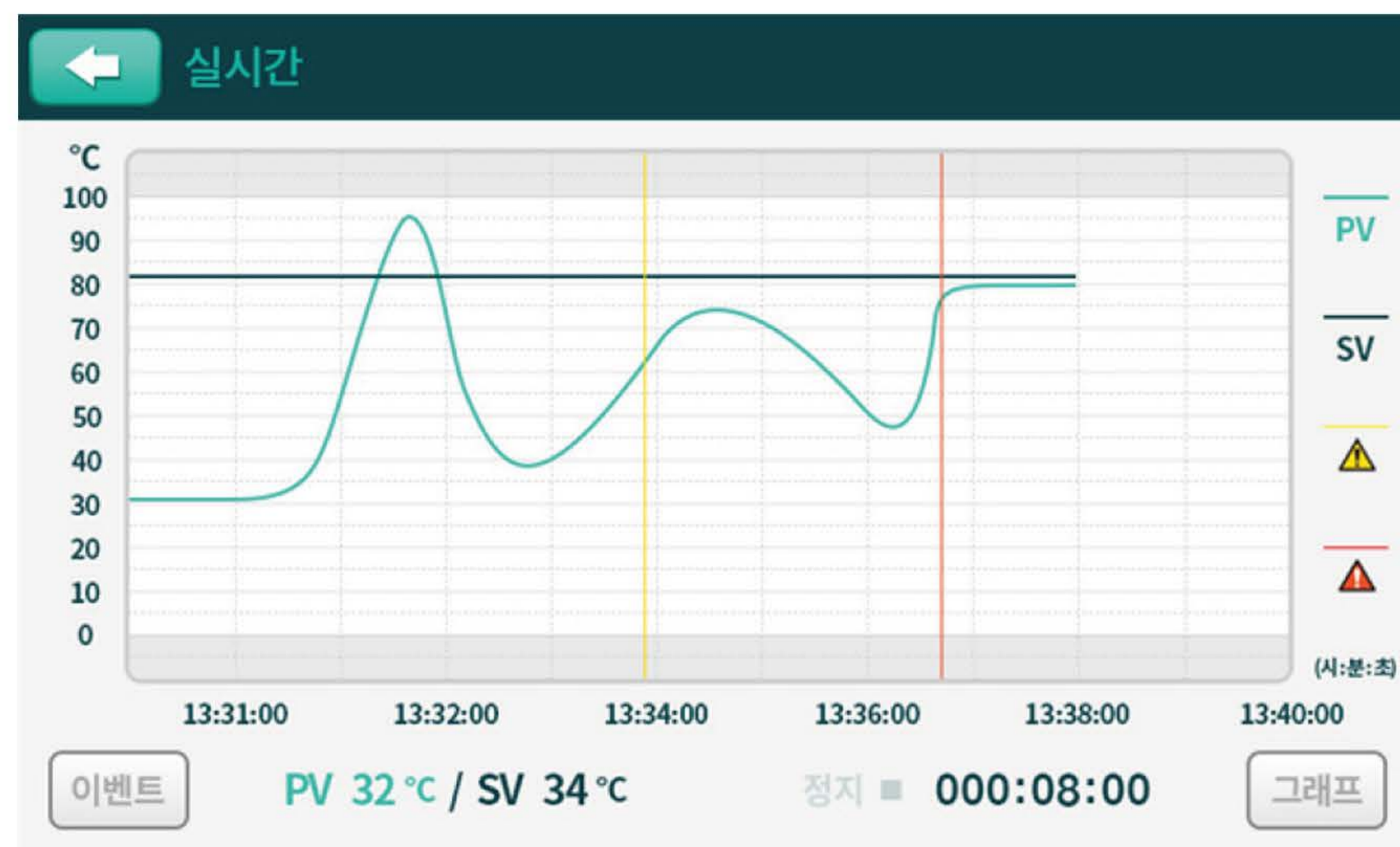

Real-time graphs that allow users to view events and alarms at a glance so that the test progress can be easily and intuitively monitored.




Test data is automatically saved(54 units) for 1 month and easy to edit and use via USB (in CSV format)

4 Main Features for New OF4 Series

04

Safe Release of Fumes

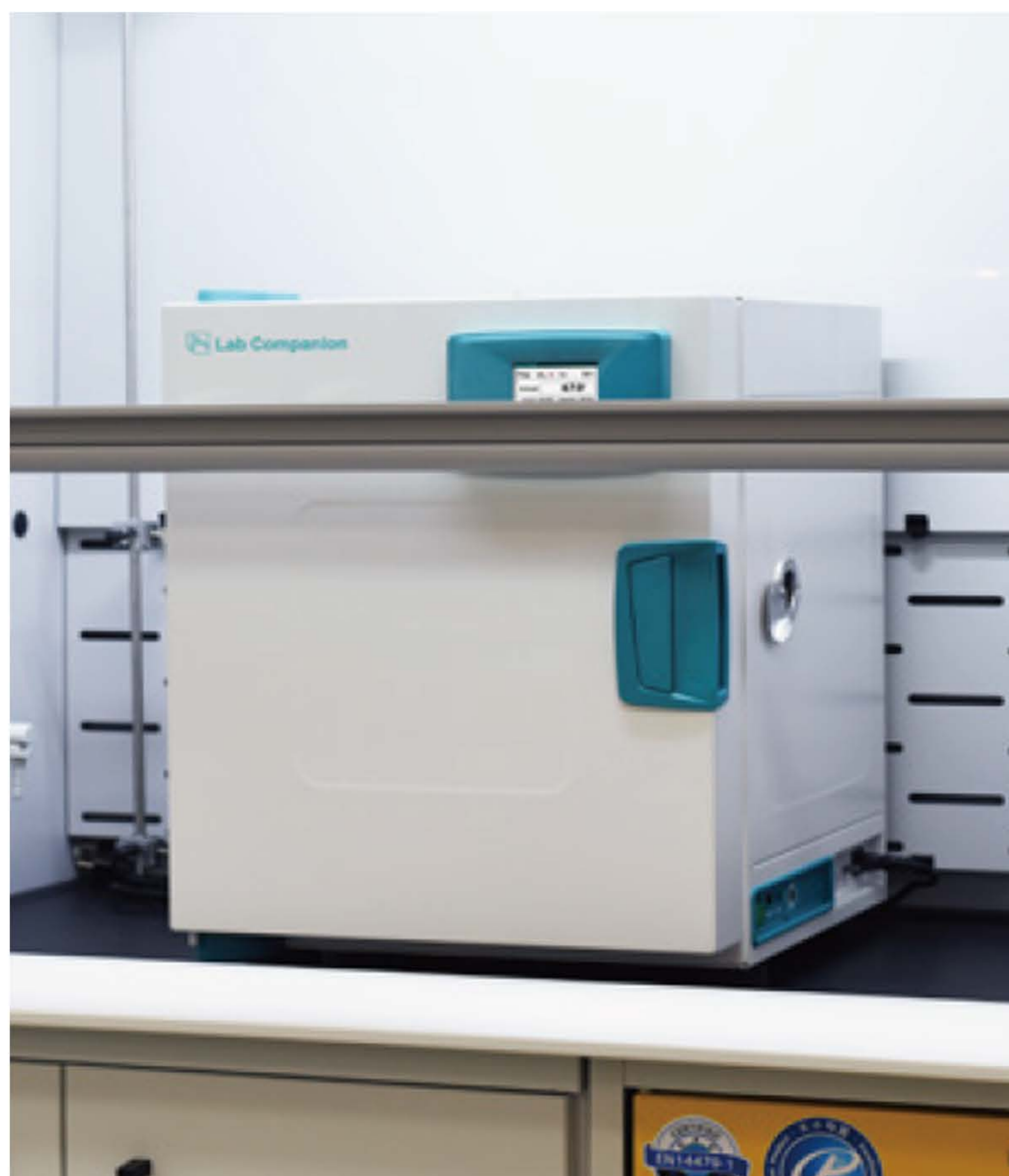

Ventilation inside a fume hood (OF4-03)



Ventilation using an arm hood

- Compact models designed to be used in a fume hood available. (OF4-03)
- Upper vent hole that can be connected to local ventilation equipment.