

Dry Cabinet

Auto desiccator



Storage of moisture-sensitive materials and samples



DC2-21L
with Shelves 4 ea (standard)



DC2-41
with Shelves 8 ea (standard)

> Automatic electronic desiccator unit

Convenient maintenance due to there being no need for a desiccant by automatically controlling the set humidity inside the storage cabinet.

> Convenient humidity setting and control

20~40%RH (Actual set-up humid.), 10~60%RH (Setting range, 1% scale)
Deviation notification and calibration for humidity & temperature allow more reliable control.

> Real-time temperature and humidity display

Displays the temperature and humidity inside the library in real time.

> Efficient two-level loading structure

It is possible to fix the feet of the product on the upper part of the product, making 2-level loading possible.
(applicable models: DC2-11(A), DC2-21S(A))

> Anti static model

Anti-Static protection via special paint the exterior
(DC2-11A/21LA/21SA/41A)



Current temperature value

Current humidity value

Humidity setting



Efficient two-level loading structure
(DC2-11, DC2-11A, DC2-21S, DC2-21SA)



Convenient internal observation

Convenient internal observation with large transparent window made of acrylic material.



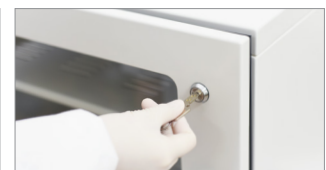
Height-adjustable shelves & optional drawers

Shelves can be adjusted in height, making it convenient loading of various samples. If you install drawers (optional), it is easy to handle storage inside.





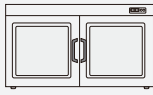
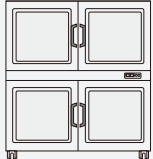


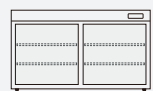
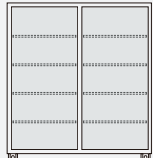
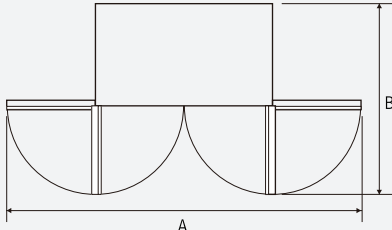
Easy to move and install

Casters which make it easy to move/fix are provided as standard for convenient moving and installation of the product. (DC2-41(A) model)



Door lock

Key lock is provided as standard to block unauthorized access.

Model		DC2-11(A) ²⁾	DC2-21L(A) ²⁾	DC2-21S(A) ²⁾	DC2-41(A) ²⁾															
Setting range ¹⁾		10~60%RH (1% scale)	10~60%RH (1% scale)	10~60%RH (1% scale)	10~60%RH (1% scale)															
Actual set-up humidity ¹⁾		20~40%RH	20~40%RH	20~40%RH	20~40%RH															
Internal volume (L / cu ft)		120 / 4.24	240 / 8.48	240 / 8.48	480 / 16.95															
Internal dimensions (W x D x H, mm / inch)		470 x 530 x 485 / 18.5 x 20.9 x 19.1	470 x 530 x 970 / 18.5 x 20.9 x 38.2	940 x 530 x 485 / 37.0 x 20.9 x 19.1	940 x 530 x 970 / 37.0 x 20.9 x 38.2															
Exterior dimensions (W x D x H, mm / inch)		555 x 715 x 650 / 21.8 x 28.1 x 25.5	555 x 715 x 1195 / 21.8 x 28.1 x 47	1115 x 715 x 650 / 43.8 x 28.1 x 25.5	1115 x 715 x 1270 / 43.8 x 28.1 x 50															
Quantity of door		1	2	2	4															
Stackable type		Yes	No	Yes	No															
Quantity of shelves (standard/max.)		2 / 4	4 / 8	4 / 8	8 / 16															
Max. load per shelf (kg / lbs)		20 / 44	20 / 44	20 / 44	20 / 44															
Weight (kg / lbs)		42 / 93	65 / 143	67 / 148	118 / 260															
Electrical requirements (230V, 50/60Hz, 1P,A)		2	2	2	2															
Cat. No.	Basic	AAHF70125K	AAHF70225K	AAHF71225K	AAHF70425K															
	Anti-electricity	AAHF70135K	AAHF70235K	AAHF71235K	AAHF70435K															
Product appearance																				
Product interior																				
		<table border="1"> <thead> <tr> <th>Model</th> <th>A (mm / inch)</th> <th>B (mm / inch)</th> </tr> </thead> <tbody> <tr> <td>DC2-11(A)</td> <td>1094 / 43.1</td> <td>1202 / 47.3</td> </tr> <tr> <td>DC2-21L(A)</td> <td>1094 / 43.1</td> <td>1202 / 47.3</td> </tr> <tr> <td>DC2-21S(A)</td> <td>2200 / 86.6</td> <td>1202 / 47.3</td> </tr> <tr> <td>DC2-41(A)</td> <td>2200 / 86.6</td> <td>1202 / 47.3</td> </tr> </tbody> </table>		Model	A (mm / inch)	B (mm / inch)	DC2-11(A)	1094 / 43.1	1202 / 47.3	DC2-21L(A)	1094 / 43.1	1202 / 47.3	DC2-21S(A)	2200 / 86.6	1202 / 47.3	DC2-41(A)	2200 / 86.6	1202 / 47.3		
Model	A (mm / inch)	B (mm / inch)																		
DC2-11(A)	1094 / 43.1	1202 / 47.3																		
DC2-21L(A)	1094 / 43.1	1202 / 47.3																		
DC2-21S(A)	2200 / 86.6	1202 / 47.3																		
DC2-41(A)	2200 / 86.6	1202 / 47.3																		

- 1) Guarantee 20 to 40% RH control within humidity setting range(10 to 60%RH) due to a deviation depending on the installation environment.
 2) DC2-A model is an anti-static product that prevents static electricity.

Lab Companion Korea

Head office & factory

Add. 153, Techno 2-ro, Yuseong-gu, Daejeon, 34025, South Korea

International sales office

Add. 10F-1005, 219, Gasandigital 1-ro, Geumcheon-gu, Seoul, 08501, South Korea

Tel. +82 2 2627 3816

Fax +82 2 3143 1824

E-mail overseas@jeiotech.com

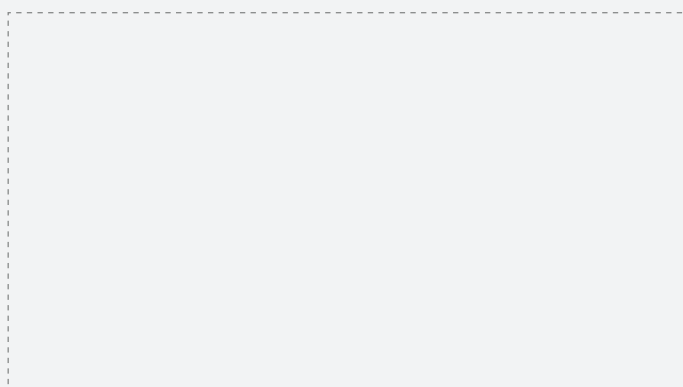
Lab Companion U.S.A. - Jeio Tech, Inc.

Add. 19 Alexander Road, Ste. 7, Billerica, MA 01821-5094, U.S.A.

Tel. +1 781 376 0700 Fax +1 781 376 0704

E-mail info@jeiotech.com

Your best business partner



www.JeioTech.com